



KLOW Blend

Results and chromatograms

Sample vial 6, injection 1 | Client report CBI-KLOW-80-091826

Laboratory Analysis: LAB-CA-02

Laboratory identity and applicable credentials are maintained by Testides and are available to regulatory authorities where required by law.


VIAL PHOTOGRAPH	SAMPLE INFORMATION	ANALYSIS	RUN DETAILS
	Peptide: KLOW Blend	Method: HPLC-UV	Run date: Sep 23, 2026
	Client ID: CBI-KLOW-80-091826	Channel: 214 nm	Acquired: 20:39:47 UTC
	Label claim: 80 mg	Integration: Laboratory	Run time: 25 min
	Injection: 10 uL, 214 nm	Identity: Lab-assigned peak	Instrument: HPLC 01

RESULTS AT A GLANCE

PROJECTED AMOUNT 74.44 mg* SUM OF COMPONENTS; SOME EXTRAPOLATED	HPLC PURITY 99.88% NAMED COMPONENTS, 214 NM*	CONTENT VS. LABEL 93.0% 74.44 MG OF 80 MG CLAIMED	IDENTITY Assigned LAB-ASSIGNED PEAK**
GHK-CU 44.10 mg* 88.2% OF 50 MG LABEL; EXTRAPOLATED	KPV 10.24 mg† 102.4% OF 10 MG LABEL	BPC-157 8.71 mg† 87.1% OF 10 MG LABEL	TB-500 11.39 mg† 113.9% OF 10 MG LABEL

CHROMATOGRAM

Every plotted reading as recorded - overlapping component peaks drawn as peak shapes fitted to the trace

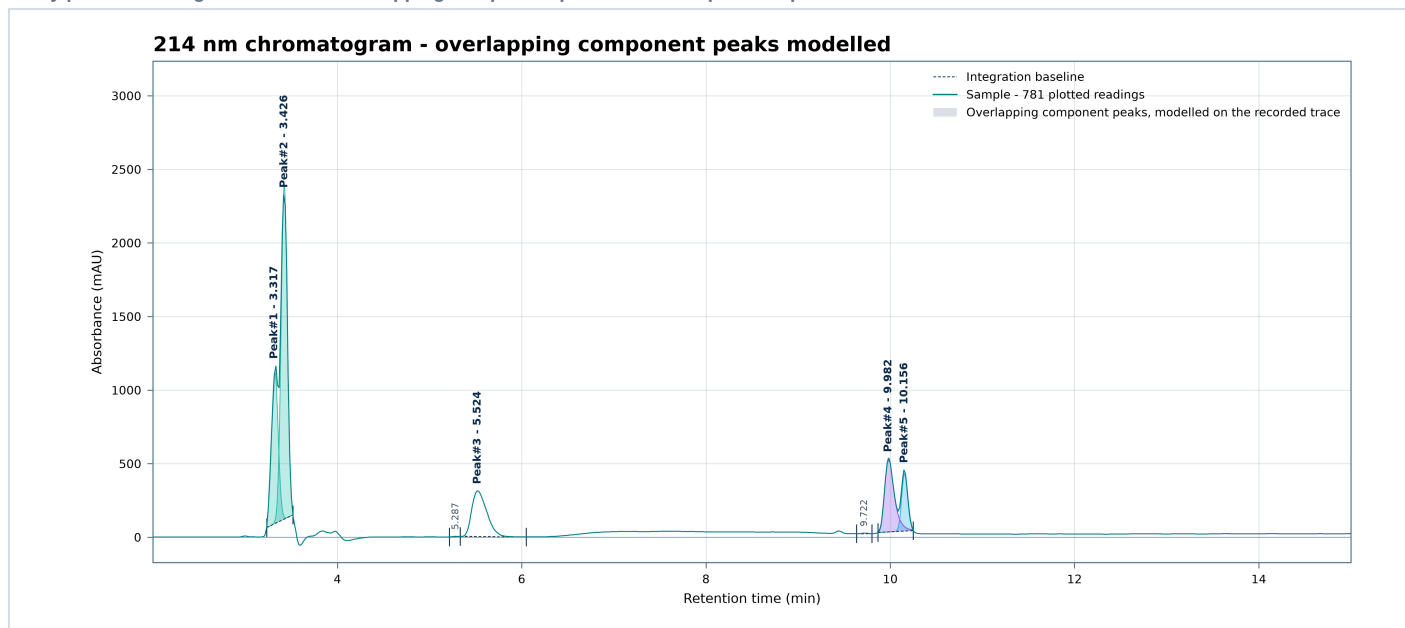

HPLC PURITY (214 NM)

Named components 99.88% Laboratory integration	Integrated peaks 7 Counted in the total	Smallest peak integrated 0.04% Of total area	Recovered chromatogram ±2.10% Area vs. laboratory table
---	--	---	--

Purity is the named component peaks' share of the total peak area the laboratory integrated at 214 nm. The smallest peak it integrated was 0.04% of the total; the laboratory does not state a reporting threshold, so peaks below the ones it integrated are not counted. Integrating the chromatogram recovered from the laboratory's report reproduces its component area to within 2.10%. This result rests on a single injection.

The chromatogram is recovered point for point from the laboratory's report: 781 plotted readings from 2.00 to 15.00 min, joined as drawn; nothing is smoothed or resampled.

* HPLC purity is the named component peaks' share of all peak area the laboratory integrated; every peak is itemised on page 3.

** Identity: the laboratory's processing method names the analyte's peak. The laboratory holds no reference standard for this analyte, so retention time is not compared.

LAB ID: D







KLOW Blend

Sample vial 6, injection 1 | Client report CBI-KLOW-80-091826

Reference comparison and calculation



CONCENTRATION CALIBRATION - 214 NM (LABORATORY VALIDATION CURVES)

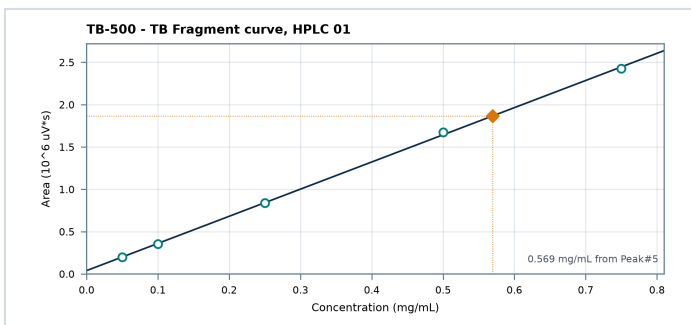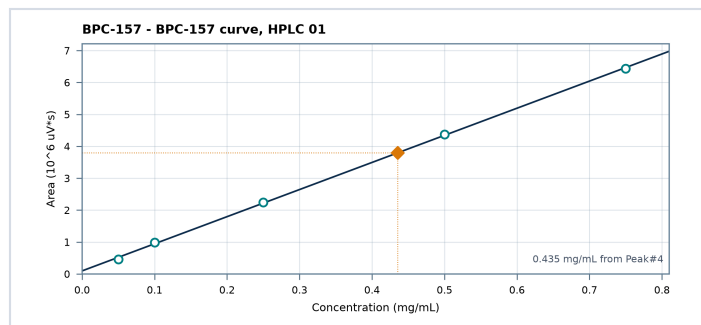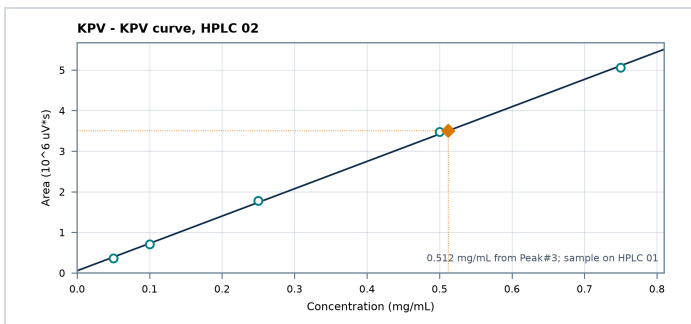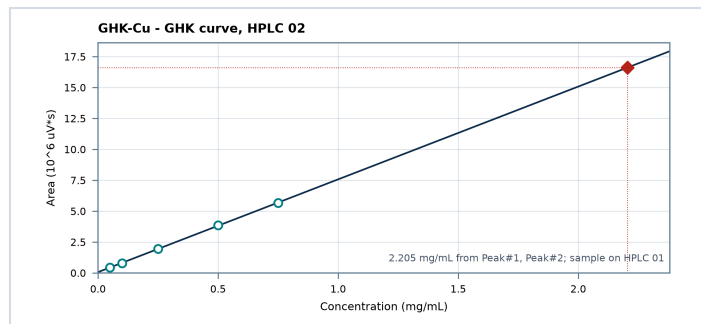

Component	Peaks	Area %	mg/mL injected	Laboratory curve	Amount	Label share	Content
GHK-Cu	Peak#1, Peak#2	54.29%	2.2051 *	GHK, HPLC 02	44.10 mg	50 mg	88.2%
KPV	Peak#3	21.27%	0.5120	KPV, HPLC 02	10.24 mg	10 mg	102.4%
BPC-157	Peak#4	15.17%	0.4355	BPC-157, HPLC 01	8.71 mg	10 mg	87.1%
TB-500	Peak#5	9.15%	0.5694	TB Fragment, HPLC 01	11.39 mg	10 mg	113.9%

† The laboratory does not report its preparation; the vial is taken as made up to 20 mL, which puts the component present in the smallest amount at the working 0.5 mg/mL. Label shares follow the blend's recorded composition: GHK-Cu 50 mg, KPV 10 mg, BPC-157 10 mg, TB-500 10 mg. Total 74.44 mg. * Outside the range of the laboratory's curve; the reading is extrapolated. Curves validated on HPLC 02 are applied to this HPLC 01 injection.

NAMED COMPONENT PEAKS

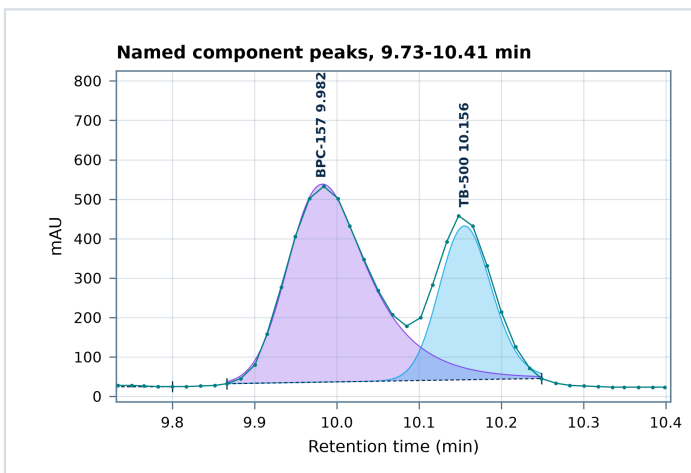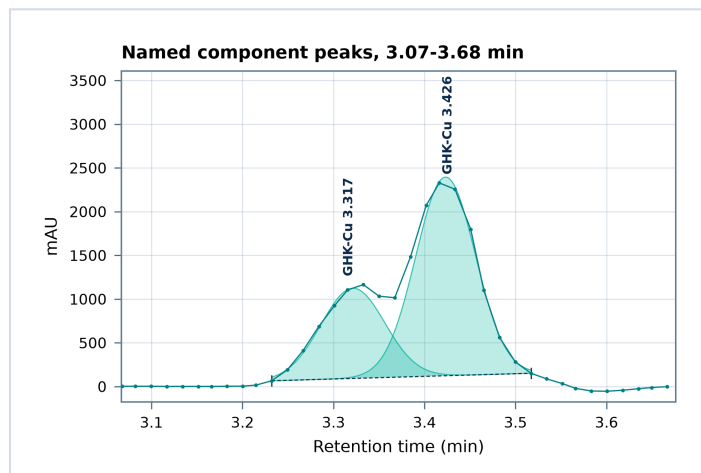

Native plotted readings, nothing smoothed. Overlapping component peaks are drawn and read as peak shapes fitted to the readings, on one baseline under the group, not the laboratory's valley-to-valley integration, which leaves the area beneath the valley in neither peak. Purity is the laboratory's. Components are assigned in the order they elute on the laboratory's method; a component can show as more than one peak.



KLOW Blend

Sample vial 6, injection 1 | Client report CBI-KLOW-80-091826

Definitions, calculation detail, and scope

HOW TO READ THIS REPORT

Start with the result on page 1. Pages 2 and 3 show the curve, comparison, and peaks behind it.

01	Recovered chromatogram The laboratory's trace, read back point for point from its report.	02	Integrated peak table Every peak the laboratory integrated, and which one is the analyte.	03	Validation curve Converts the peak area to a concentration for this analyte.
-----------	---	-----------	---	-----------	--

WHAT THE RECOVERED CHROMATOGRAM ADDS

How the chromatogram on this report was obtained	
How it was recovered The laboratory's data system draws the chromatogram in its report as vector data. Every vertex it drew was read back - 781 readings - and both axes were calibrated from the report's own tick marks. Integrating the recovered trace over the laboratory's baselines and peak marks reproduces the area in its peak table to within 2.10%.	What it can and cannot show It shows every peak, shoulder and baseline the laboratory integrated, at the resolution the laboratory printed. It carries a single wavelength and no UV spectrum, so it cannot reveal a co-eluting compound the way diode-array data can.
Scope: the recovered chromatogram reproduces the laboratory's report. It does not re-measure the sample.	

HOW THE PURITY RESULT WAS MEASURED

214 nm area percent, from the laboratory's own integration							
The laboratory integrated each peak to its own baseline and named the analyte's peak in its processing method. Purity is that peak's share of the total integrated area; the shares below are recomputed from the areas, and each peak's area is checked against the recovered chromatogram.							
Purity (%) = assigned peak area / total integrated peak area x 100							
Peak	RT (min)	Area (uV*s)	Area %	Tailing	Plates	Role	Recovered
1	3.317	1,736,401	10.54	0.80	13,653	Peak#1 - GHK-Cu	1,724,208 (-0.70%)
2	3.426	7,204,519	43.74	1.18	12,464	Peak#2 - GHK-Cu	7,053,038 (-2.10%)
3	5.287	6,952	0.04	0.85	41,311	impurity	-
4	5.524	3,502,906	21.27	1.39	5,213	Peak#3 - KPV	3,517,575 (+0.42%)
5	9.722	12,595	0.08	0.99	105,179	impurity	-
6	9.982	2,499,191	15.17	0.99	57,350	Peak#4 - BPC-157	2,519,230 (+0.80%)
7	10.156	1,507,534	9.15	1.14	98,842	Peak#5 - TB-500	1,496,572 (-0.73%)
Integrated peaks 7 smallest 0.04% of total injections 1							
Area percent assumes every peak absorbs at 214 nm like the main peak; the analyte's peak is the one the laboratory's method names. Recovered areas are compared for the analyte's peaks of 1% or more.							

WHAT THIS REPORT CAN AND CANNOT TELL YOU

The report explains the measurement; it does not decide how a material should be used	
Observed evidence The laboratory's chromatogram, its integration, and its validation curve.	Identity and composition No definitive identity or absolute purity unless separately supported.
Calculated results Area-percent purity, concentration, projected amount, and calibration trace.	Safety and performance No conclusion on safety, sterility, efficacy, fitness for purpose, or suitability.
Method context Laboratory peak assignment; retention against its reference standard where methods match.	Regulatory status No regulatory compliance finding, approval, or authorization to manufacture, sell, or use.
This analytical record does not authorize or endorse manufacture, sale, marketing, distribution, or use of the material.	